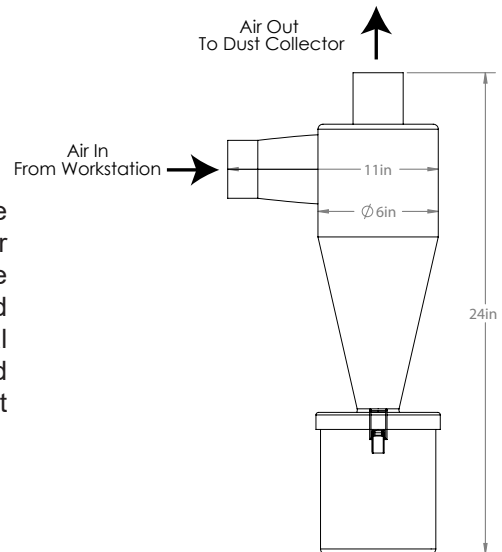

THE 3 " ACCUMULATOR

- 1 - Accumulator
- 1 - Mounting Bracket Kit
- 1 - Large Circular Clamp
- 2 - #8 Screws

INSTALLATION AND MOUNTING :

The Accumulator comes with a mounting bracket that may be fastened to a wall by using the two #8 wood screws provided or other wall fasteners to attach the bracket in place before assembling the Accumulator. To fasten to a wall, loosen the large band clamp and remove from the Accumulator. Fasten to the wall or other vertical surface before sliding the Accumulator back into the clamp and tightening. Leave enough room (1-2 inches) under the cup so that it may be emptied.



HOSE CONNECTION:

The Accumulator is installed between the workstation and the vacuum system. Connect the hose from the dust collector to the outlet on top of the Accumulator. The hose collecting the debris should connect to the inlet tube, horizontal arm, of the Accumulator.

The inlet and outlet tubes of the Accumulator are designed to accommodate 2- 1/2" shop-vac hose or standard 3" vacuum hose. The 2-1/2" hose end connectors slip inside the inlet and outlet tubes and seal with the foam rubbergaskets supplied. Insert the hose end by gently pushing and twisting it into place. For standard 3" vacuum hose, remove the foam rubber gasket inside the tube, then slip the hose over the inlet or outlet tube and clamp with a hose clamp.

OPERATION:

Make sure the collection cup is snapped into place on the bottom of the unit. The metal toggle latches hold the cup in place. To remove the cup, make sure the vacuum system is off. Gently tap the sides of the Accumulator to shake down dust into the cup. Hold the collection cup with one hand on each side directly below the toggle latches. With your index fingers, flip the toggles upward then swing the wire bales outward. Lower the cup and empty the contents. Replace the cup onto the bottom rim of the Accumulator, swing the wire bales over the bottom rim of the Accumulator, and snap the toggle latches down to secure the cup in place.

Turn the vacuum system on and use as you normally would. To maintain excellent suction, it is best to check the cup daily and empty it when it is half full.