

Vanguard Gold Mobile

Extraoral Chairside Suction




Part
Number:
10337



HEPA filtration of viruses, chairside biological aerosols, and moisture droplets.



Long 6ft flexible hose can be adjusted and stays in place.

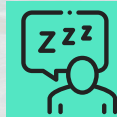

Super Quiet
Operation

Adjustable
Speed



Long lasting
powerful motor.

Easy to move
and operate.



Optional on
& off foot
pedal.



Product Specs



120VAC 15A/5A
(220V/240V) Opt.



L 14"xW 10.5"xH 14.5"



WARRANTY



39 LBS

- Internal cartridge HEPA filter
 - Filters 99.97% of particles down to 0.3 micron
- Internal particle filter bag
 - Filters down to 1 Micron
- Rolling caster wheels allows easy movement
- Flex hose links can be added/removed & disinfected
- Optional exhaust mounted Odor filter
 - Can absorb mercury vapors from amalgam filling removal
- Replace bag filter as needed as it will be replaced most frequently
 - check 1 time per 3 months.
- Replace HEPA filter as needed based on appearance
 - check 1 time per 6 months